INN
911

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IMPORTANT SAFETY INSTRUCTIONS:

- Read the operating instructions carefully before use and save this user manual for future references.
- This device is intended for household use only.
- This product can be used by children over 8 years old and by persons with reduced physical, mental or sensory abilities, or lack of experience and knowledge, only if under supervision or if they have received instructions regarding the safe use of the appliance and understand the associated dangers.
- Make sure children do not play with the appliance. Cleaning and maintenance operation should not be performed by children without supervision.
- Never immerse the appliance or its base in water or other liquids.
- When the charger is connected to the mains do not handle the appliance with wet hands or feet.
- Only use the accessories recommended by the manufacturer.
- Do not use the appliance if the main unit or the brush head are damaged.
- If more than one person uses the toothbrush, we recommend for hygienic reason that each person use their own brush head.
- Do not pull the power cord to remove the appliance.

- Do not wrap the power cord around the appliance. Do not fold or twist the power cord.
- Unplug the appliance in case of a power failure and before proceed with any cleaning and maintenance operations.
- Do not dismantle or try to repair the appliance. The repairs should be performed by the authorized after sales service or qualified personnel.
- This product is not a toy.
- Keep this product away from heating sources.
- The power cord can not be replaced. If the cable is damaged, the appliance must be replaced.
- This product is not intended for professional use in dental offices, only use the appliance to clean teeth and gum. In case of pain, please contact your physician.
- The battery included in the brush is not removable by the user and will be the recovery and treatment center to extract it before disposal.

TECHNICAL SPECIFICATIONS

Battery Type: Lithium battery 3.7V/ 400mAh

Rated Voltage: 220-240V ~ 50/60Hz

Power: 2W

Waterproof: IPX7


Spare toothbrush, model INN-90801

Charging Base

LED indicator for charge level

PRODUCT DESCRIPTION

1. Holder for brush head
2. ON/OFF and SPEED button
3. LED indicator for charge level
4. Bristles
5. Brush Head
6. Charging Base
7. Power Cord.

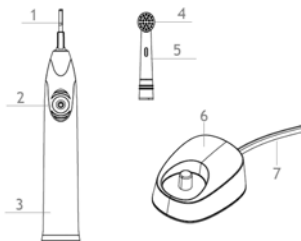HOW TO USE

Charging the battery

1. Insert the plug in a power outlet.
2. Place the toothbrush on the charging base. The LED indicator becomes red.
3. It takes at least 7 hours to complete the charging cycle. The LED indicator becomes green when the battery is fully charged.
4. When the battery is fully charged, you can remove the plug from the power outlet.
5. Charge the battery when you note a significant decrease in speed, or when the brush head is completely stopped and the LED indicator becomes red.

NOTE:

- . Turn the toothbrush off during the charging base.
- . Charge the battery before using the toothbrush for the first time.
- . The charger gets warm during the charging. This is quite normal.
- . The battery must be charged with the charger recommended by the manufacturer.

HOW TO REPLACE THE BRUSH HEAD

1. Remove the head.
2. Insert the head directly into the holder as shown in the following figure.

USE OF THE PRODUCT

1. Put the required amount of toothpaste on the toothbrush.
2. Put the brush head in your mouth and rub the brush against your front teeth to press the toothpaste between the bristles.
3. Press the power switch to turn the toothbrush on and start brushing at the back of your mouth.
4. Brush over all sides of your teeth, between your gums and between your teeth. 2-3 minutes of brushing are usually sufficient.
5. After brushing your teeth, rinse the brush head under running water.

NOTE: Do not turn the toothbrush on before putting it in your mouth and switch it off before taking it out to prevent splashing.

TIMER: Normal tooth cleaning requires at least 2-3 minutes. The appliance is equipped with a useful timer. After 2 minutes the head stops for a moment, and then continues to rotate.

CLEANING AND MAINTENANCE


After use, rinse the head under running water, remove excess water and store it in an upright position.

Do not immerse the toothbrush or its base in water. Remove the brush head and clean the main unit with a damp cloth.

It is recommended to replace the brush head every 3 months and every time you notice an excessive consumption of the bristles. Use only the heads suitable for this model.

After use store the product in a dry place away from heat sources.

RATING LABEL:

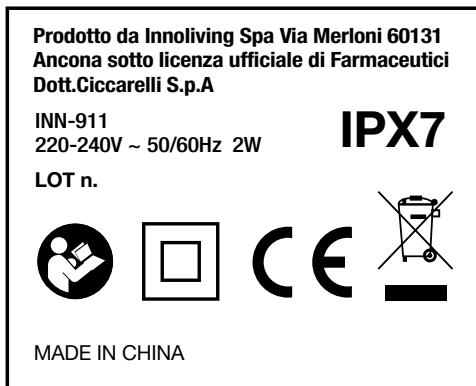Read carefully the instruction manual.

IPX7

Degree of protection from contact with water.

Class II insulated electrical appliance

This product complies with all the applicable EU Directives.

INFORMATION TO USERS according to Legislative Decree No. 49 of March 14, 2014 “ 2012/19/UE Implementation of the Directive on Waste Electrical and Electronic Equipment (WEEE) “The crossed-out dustbin symbol indicates that the product at the end of its life must be collected separately from other waste.

The user should, therefore, take the equipment with the essential components at the end of its useful life to the separate collection center of electronic and electrical waste, or return it to the retailer when purchasing new equipment of equivalent type, in ratio of one to one, or one to zero for devices with larger side less than 25 CM. The separate collection for the decommissioned equipment for recycling, treatment and environmentally compatible disposal contributes to avoid possible negative effects on the environment and human health and promotes recycling of the materials. Improper disposal of the product by the user entails the application of administrative sanctions according to Legislative Decree No. 49 of 14 March 2014

The batteries used in this device must be disposed of in the special bins at the end of their life. Please inform yourself about the local rules on separate collection of batteries.

The correct disposal of batteries helps preventing potentially negative consequences on the environment and human health.

innoliving
— — *easy tech, easy life*

Prodotto Da Innoliving Spa
Via Merloni 9 60131 Ancona
Tel: 071 2133550
www.innoliving.it

MADE IN CHINA