INN  
908

Testina  
compatibile  
**ORAL-B®**

**ORAL-B®**  
compatible  
brush head

Glava  
kompatibilna s:  
**ORAL-B®**

Zdržljiva  
glava ščetke  
**ORAL-B®**

SPAZZOLINO ELETTRICO / TOOTHBRUSH MODEL / ELEKTRIČNA ČETKICA ZA ZUBE / ELEKTRIČNA ZOBNA ŠČETKA

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#### EN

*Pasta del Capitano toothbrush model INN-908 is endowed with a circular head able to surround and wrap each single tooth and help to remove dental plaque more effectively.*

#### TECHNICAL SPECIFICATIONS

- Rotating Speed: 6000 times/minute
- Rotating angle: 15 degree
- Filament: soft, end-rounded bristle
- Replaceable head
- Battery Type: 2 x AA battery included

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2 x 1.5V AA carbon zinc batteries (included)

Modello / Model: R6P - Fabbricata da / Manufactured by: Zhejiang Mustang Battery Co., Ltd. NO.818 Rongji Road, Luotuo, Ningbo, 315202 Zhejiang, P.R. China

Importata da / Imported by: Innoliving Spa Via Merloni 2/B - 60131 Ancona, Italia

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*Read the following instructions and warnings before use and keep them as reference for future uses.*

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#### WARNINGS

*Use the product for the intended purposes only, that is to say as electric toothbrush. This product is not recommended for use by children under age 3. This product has been conceived for HOUSEHOLD use only, do not use the device for professional uses. Keep it far away from heating sources. Ask for your dentist's advice on the way and frequency of use. Use it only for cleaning tooth and gums. In case of pain or bleeding, please contact your doctor.*

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#### BATTERY REPLACEMENT


Open the battery room, remove and replace the batteries according to the indicated polarities. The toothbrush works with 2 AA batteries. Remove the batteries from the use in case of long no-use.

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#### HOW TO USE

Turn ON/OFF the unit, pressing the ON/OFF button. After use, rinse the brush head under tap water, eliminating the excess of water and store it vertically. Do not immerse the handle of toothbrush into water, for the cleaning, remove the brush head and clean the handle with a damp cloth. It is recommended to replace the brush head every 3 months and every time an excessive consumption of the bristles is noted. Use only the spare brush head models INN-90801NEW, INN-90802, INN-90803, INN-90804, INN-90805, INN-90806 produced for this item.

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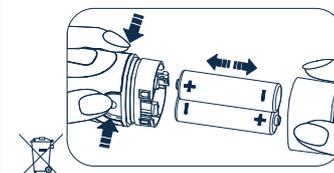#### INFORMATION TO USERS

according to Legislative Decree No. 49 of March 14, 2014<sup>9</sup> 2012/19/UE Implementation of the Directive on Waste Electrical and Electronic Equipment (WEEE)<sup>10</sup>. The crossed-out dustbin symbol indicates that the product at the end of its life must be collected separately from other waste. The user should, therefore, take the equipment with the essential components at the end of its useful life to the separate collection center of electronic and electrical waste, or return it to the retailer when purchasing new equipment of equivalent type, in ratio of one to one, or one to zero for devices with larger side less than 25 CM. The separate collection for the decommissioned equipment for recycling, treatment and environmentally compatible disposal contributes to avoid possible negative effects on the environment and human health and promotes recycling of the materials. Improper disposal of the product by the user entails the application of administrative sanctions according to Legislative Decree No. 49 of 14 March 2014

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The batteries used in this device must be disposed of in the special bins at the end of their life. Please inform yourself about the local rules on separate collection of batteries. The correct disposal of batteries helps preventing potentially negative consequences on the environment.